# Aerogen® GO

## Vibrating Mesh Nebulizer

### Instruction Manual

- EN English
- DE Deutsch
- FR Français
- ES Español
- NO Norsk
- SV Svenska
- AR عربي

Actual size

**1 English**

**17 Deutsch**

**33 Français**

**49 Español**

**65 Norsk**

**81 Svenska**

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English

---




# **Aerogen<sup>®</sup> Go Vibrating Mesh Nebulizer**

Instruction Manual

# Short Guide

## Assemble


For complete assembly of all components, please refer to **pages 6-9** of the instruction manual.
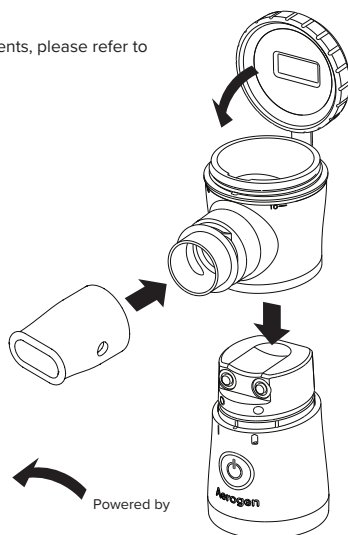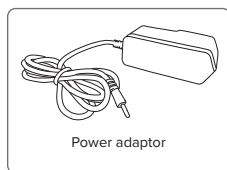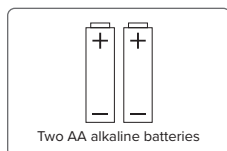
## Inhale

For more details, please refer to **page 9** of the instruction manual.

Breathe through the mouthpiece or the mask. Sit in an upright position and hold the Aerogen® Go vertically.


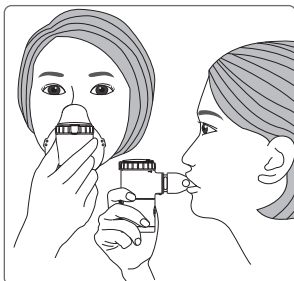
## Cleaning and disinfection

For more details, please refer to **page 10** of the instruction manual.



### • Cleaning after every operation.

- Warm filtered water with mild liquid dishwashing soap
- Rinse in filtered water
- Air dry completely on a clean, lint-free cloth.
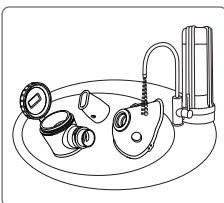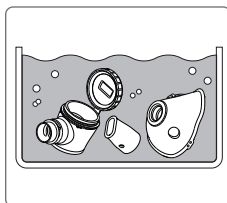
### • Daily disinfection

- Standard commercial disinfectant for baby bottles or soak in 70%-75% ethyl alcohol for 1 minute
- Air dry completely on a clean, lint-free cloth
- Wipe the main unit with a clean cloth moistened with 70%-75% ethyl alcohol.
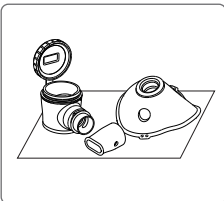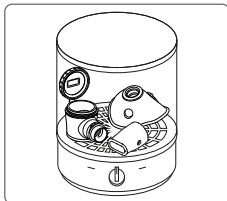
# Table of Contents

|                                                |    |
|------------------------------------------------|----|
| Important Information .....                    | 4  |
| The Aerogen® Go Vibrating Mesh Nebulizer ..... | 5  |
| Names and Functions of the Parts .....         | 6  |
| How to Set-up Your Aerogen® Go Nebuliser ..... | 7  |
| How to Use Your Aerogen® Go Nebuliser .....    | 8  |
| Cleaning .....                                 | 10 |
| Disinfection .....                             | 10 |
| Storage and Transportations .....              | 11 |
| Troubleshooting .....                          | 11 |
| Disposal .....                                 | 12 |
| Spare Parts and Accessories .....              | 13 |
| Product Specifications .....                   | 13 |
| Technical Data .....                           | 13 |
| Electromagnetic Compatibility .....            | 14 |
| Symbols and Machine Markings .....             | 16 |
| Warranty Information .....                     | 16 |

The Aerogen® Go Vibrating Mesh Nebulizer incorporates our leading-edge Micro Synchronised Delivery Technology (MSDT) to empower people with respiratory conditions, such as asthma, chronic bronchitis, cystic fibrosis, emphysema, COPD (chronic obstructive pulmonary disease), and PAH (pulmonary arterial hypertension), to control the respiratory health and simplify the treatment. Our iconic design is distinct from currently available products. Aerogen® Go is efficient, water and dust resistant, lightweight, and quiet.

# Important Information

## 1. Intended Use

- The Vibrating Mesh Nebulizer is an inhaled delivery system that aerosolises approved medicines for respiratory treatment.

## 2. Intended Users

- The nebulising system can be used by patients over the age of 6 months who require inhalation treatment and can spontaneously breathe. Children and anyone who require assistance or supervision must be observed constantly by an adult, caregiver, or medical professional when using the device for inhalation therapy.

## 3. General

- To ensure safe and correct use of the Aerogen® Go, please read this instruction manual carefully and in its entirety before using.
- Keep this manual in a safe place for future reference.
- The Aerogen® Go Vibrating Mesh Nebulizer is designed for a single patient. Do not share the Aerogen® Go with other users.

## 4. Information About the Instructions for Use

- The information contained in these instructions may be modified in accordance with legal requirements. For the most up-to-date product information, please contact [info@aerogen.com](mailto:info@aerogen.com)

## 5. Safety Instructions and Warnings

Please read the following warnings carefully and completely. Failure to observe these could result in minor to moderate personal injury, or may impair the effectiveness of your device.

- Do not share your Aerogen® Go with any other users. It is designed for single patient use only.
- Carefully and thoroughly clean and disinfect your Aerogen® Go before using it for the first time. Your Aerogen® Go must be cleaned after each use and disinfected daily after the last use of the day and cleaning procedure. Instructions for proper hygiene are printed in page 10 of this manual.
- Do not allow the main unit or the power adaptor to come in contact with water or any other liquid. If liquid touches either the main unit or the power adaptor, discontinue use immediately and wipe the liquid away.
- In case of electric shock, do not attempt to interact with the power adaptor with wet or damp hands.
- Do not place your Aerogen® Go where infants, babies, and children can access it. Keep out of reach of children at all times. Infants and children may only use Aerogen® Go when under constant adult supervision.
- Do not allow the use of Aerogen® Go on unconscious patients who cannot breathe spontaneously.
- In order to reduce the risk of damage to medicine, or infection, do not fill the medication reservoir when it is not in use. Only add medicine when you are preparing to begin inhalation therapy.
- Do not fill the medication reservoir with water, essential oils, or any non-medicinal liquids or substances. Only use medicines that have been prescribed to you by your healthcare professional. Failure to adhere to this instruction may cause health problems and could damage the nebulising system, voiding your warranty.
- Please read the following product notices carefully when using your Aerogen® Go. Failure to adhere to these could result in damage to your device.
- Do not soak your Aerogen® Go main unit or power adaptor or allow them to become submerged in water or any other liquid.
- When your inhalation therapy session is finished, switch off your Aerogen® Go immediately. Failure to do so may damage the nebulising system or other components.
- The mesh mechanism is highly sensitive to damage. Do not touch it at any time. Doing so may cause damage and may void your warranty.
- Do not drop your Aerogen® Go or cause it to impact any surface, whether it is in use or not. Doing so may cause damage and may void your warranty.
- Do not mix multiple brands or types of batteries. Do not mix old and new batteries. Replace all batteries at the same time.
- Do not remove any screws from your Aerogen® Go.

## Important Information


- Your Aerogen® Go does not have a battery charging function.
- Your Aerogen® Go will not function properly if the medication reservoir is not properly attached when the device is on.
- Do not touch the functional membrane attached to the medication reservoir cap.
- Separate the medication reservoir from the main unit before attempting the cleaning and disinfection procedure.
- Do not poke, break, or disassemble any Aerogen® Go components.

### 6. Using Aerogen® Go with Babies Over 6 Months, Children, and Persons Who Require Assistance

- Babies, children, and anyone who require assistance or supervision must be observed constantly by an adult or medical professional when using the Aerogen® Go for inhalation therapy. This helps ensure safe and effective treatment. Two masks are included in the package with your Aerogen® Go. The paediatric mask is designed for toddlers that are unable to use the standard mouthpiece. The adult mask is designed for patients who have difficulty using the standard mouthpiece.
- Always keep the Aerogen® Go and all associated parts out of the reach of infants and small children.

## The Aerogen® Go Vibrating Mesh Nebulizer


Please examine the contents of your package and compare the included components to the pictures below. If any item is missing, please contact Aerogen Customer Service immediately – [info@aerogen.com](mailto:info@aerogen.com) or +353 91 540 400.

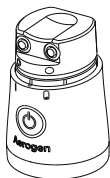
Main Unit


Medication Reservoir

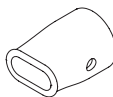Mouthpiece

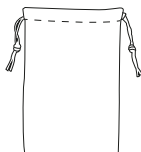
Storage Bag

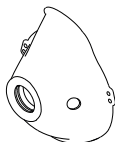Adult Mask

Paediatric Mask

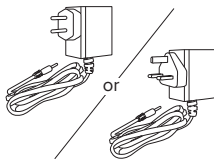
Power Adaptor  
Type C / Type G  
(optional)

Alkaline  
Batteries  
(2xAA/LR6)

# Names and Functions of the Parts

## Components

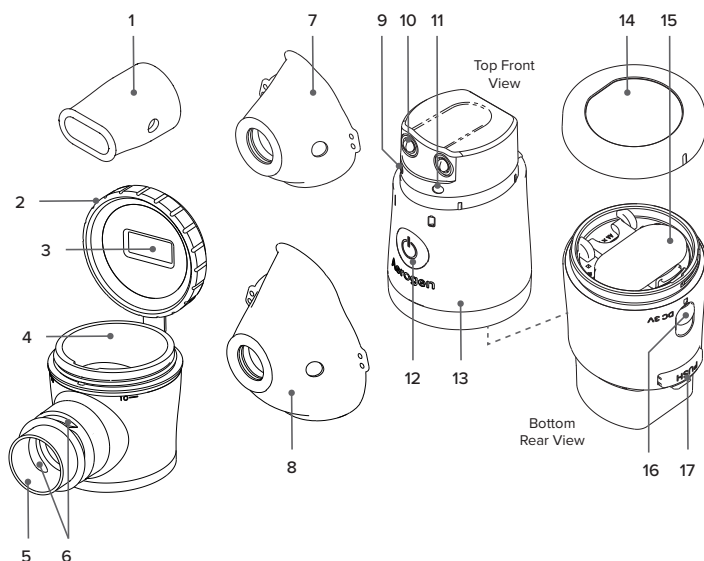1. Mouthpiece
2. Medication Reservoir Cap
3. Functional Membrane
4. Medication Reservoir
5. Aerosol Port  
(Mesh Mechanism within)
6. Vent
7. Paediatric Mask
8. Adult Mask
9. Power Indicator

10. Electrodes
11. Low Battery Indicator
12. On/Off Button
13. Main Unit
14. Outer Battery Cover
15. Inner Battery Cover
16. Power Adaptor Port/Cover
17. Medication Reservoir  
Release Button

## Operating Life

| Component                      | Expected lifetime |
|--------------------------------|-------------------|
| Main unit                      | 3 years           |
| Medication reservoir           | 1 year            |
| Mouthpiece                     | 1 year            |
| Paediatric mask and adult mask | 1 year            |


- The Aerogen® Go nebuliser is designed to be used for 10 minutes per session, 3 sessions per day. The life cycle of the components may be shortened if the nebuliser is used for longer durations.

# How to Set-up Your Aerogen® Go Nebuliser

The Aerogen® Go can be operated with batteries or power adaptor.


## Inserting Alkaline Batteries

1. Twist off the outer battery cover and open the inner battery cover.
2. Insert two AA alkaline batteries by aligning with indicated polarities.  
 ⚠ If batteries are inserted with incorrect polarities, Aerogen® Go will not power on.
3. Press down the inner battery cover until it clicks. Twist the outer battery cover to seal it.
4. When the low battery indicator flashes, battery power is insufficient. Replace the batteries with new ones as soon as possible. When the low battery indicator is lit constantly, battery power is too low to power on the device. Replace the batteries with new ones immediately.  
 ⚠ Battery life for operation could vary among different brands or types of batteries. Keeping spare batteries available is recommended.
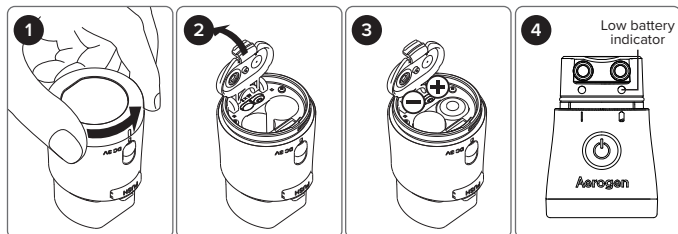
## Connecting to the Power Adaptor


1. Open the power adaptor port cover in the back of the main unit, and connect the power adaptor to this port.
2. Plug in the power adaptor to a suitable power socket.  
 ⚠ The Aerogen® Go does not have a battery charging function.  
 ⚠ Before using the power adaptor, remove the batteries from the main unit.

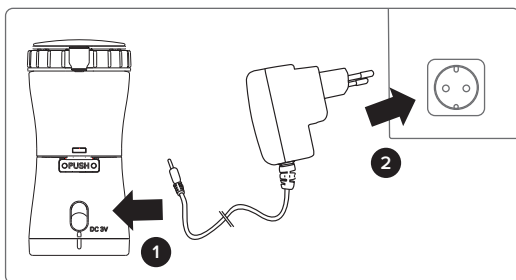
# How to Use Your Aerogen® Go Nebuliser

## Adding Medication

1. Align the base of the medication reservoir with the top of the main unit and press the two pieces together until the medication reservoir release button clicks.
  - 1.1. Before attaching the medication reservoir onto the main unit, place the main unit with the front side face to you and visually check if the silicone ring is in position.

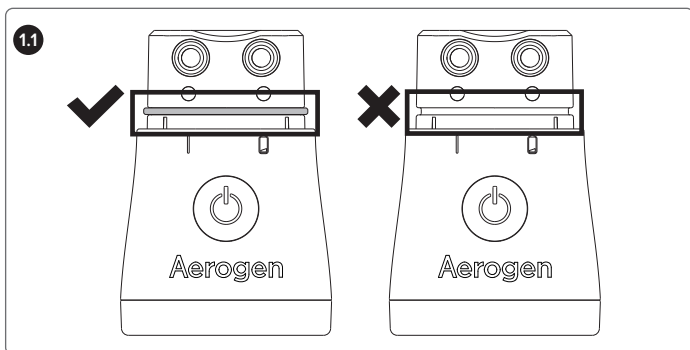
- 1.2. After attaching the medication reservoir onto the main unit, check the side face of the device. To ensure the medication reservoir is securely attached on the main unit, there shall be no significant or declined gap between the medication reservoir and the main unit.

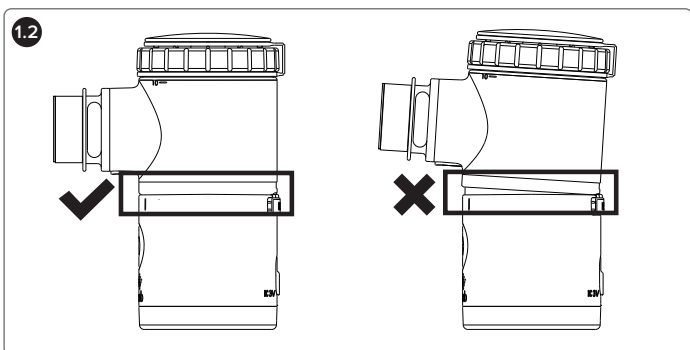
# How to Use Your Aerogen® Go Nebuliser

2. Twist the medication reservoir cap counterclockwise to open.
3. Pour the recommended amount of medicine into the open medication reservoir. The amount of medicine will vary depending on your prescribed dosage. The minimum fill volume is 1ml. Do not fill the medication reservoir with more than 10ml of medicine, as indicated by the top line on the side of the medication reservoir.
4. Twist the medication reservoir cap clockwise to close until it does not turn any more. The Aerogen® Go is now ready for use.
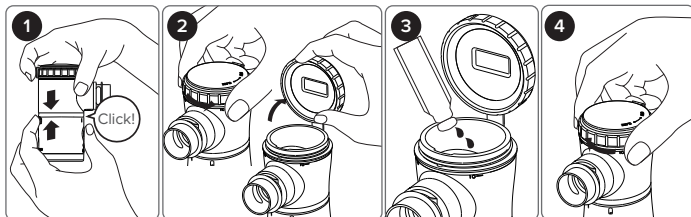
## Inhalation Therapy Session




1. Connect either the mouthpiece or mask to the medication reservoir.
2. Press the On/Off button to start the nebulisation process. When the green light is on, nebulisation has begun.
  - ⚠ If the medication reservoir is not properly attached to the main unit, the device will not turn on.
3. Sit in an upright position and hold the Aerogen® Go vertically. Breathe in. Breathe normally in order to achieve the best result.
  - ⚠ During inhalation process, please do not cover the vent holes on the mouthpiece.
4. End of the inhalation session. Press the On/Off button to turn off the device. When the medicine is running out, the aerosol output will decrease. Please switch off the power at this time.
  - ⚠ To avoid malfunction of the nebuliser, switch off the device when the nebulisation session ends or when the liquid medication in the medication reservoir is running out.
  - ⚠ For safety reasons, the device will switch off automatically after 10 minutes of nebulisation. If additional treatment is necessary, press the On/Off button to start the next nebulisation session.
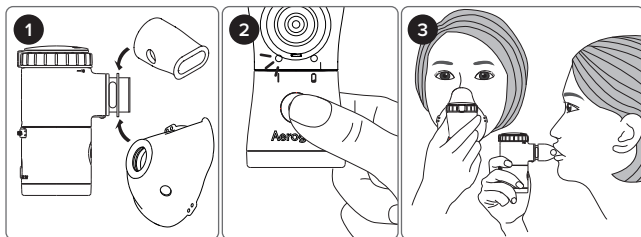
# Cleaning

## Cleaning After Every Operation

1. Remove the mouthpiece or mask from your Aerogen® Go.
2. Separate the medication reservoir from the main unit and discard residue liquid in the medication reservoir.
3. Place the medication reservoir, mouthpiece or mask in warm filtered water with mild liquid dishwashing soap and shake a couple of times.
4. Wash components mentioned above thoroughly in filtered water to remove soap residue.
5. If any components mentioned above are still stained after cleaning, repeat step 3 and step 4 until they are clean.
  - △ Do not touch the mesh with any objects or fingers. Do not soak the main unit or the power adaptor in any liquid. Only clean those components by wiping them off with clean paper towel or cloth.
6. Remove excess water by shaking components carefully and let them air dry completely on a clean, lint-free cloth.
  - △ Do not place your Aerogen® Go in an oven, a microwave or a dishwasher.
7. Place your Aerogen® Go components in the storage bag or other dry and clean containers.
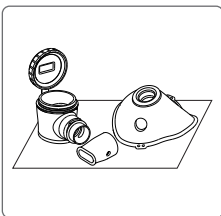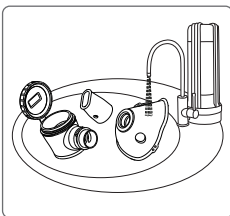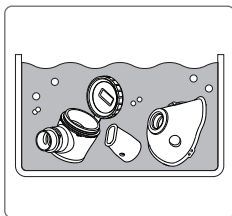
# Disinfection

## Daily Disinfection

### 1. Disinfection of the medication reservoir, mouthpiece and paediatric/adult mask

At the end of the day, after your Aerogen® Go has been cleaned as the cleaning procedure above, use one of the recommended disinfection instructions as follows.



- **Using a standard commercial disinfectant for baby bottles**

Follow the instruction provided by the manufacturer of the commercial disinfectant for an effective disinfection. Do not use an oven, a microwave or a dishwasher to disinfect your Aerogen® Go.

- **Using alcohol**

Immerse your medication reservoir, mouthpiece and paediatric/adult mask in 70%-75% ethyl alcohol for approx. 1 minute.

Remove excess liquid by shaking components carefully and let them air dry completely on a clean, lint-free cloth. Place your Aerogen® Go components in the storage bag or other dry and clean containers.
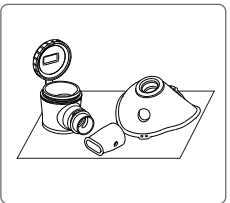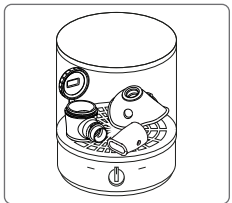
### 2. Disinfection of the Main Unit

Wipe the main unit with a clean cloth moistened with 70%-75% ethyl alcohol.

# Storage and Transportations

Remove batteries and separate medication reservoir from main unit, if your Aerogen® Go is expected to be in storage for a prolonged period. After cleaning and disinfection, place your individual components of Aerogen® Go in the storage bag or keep in a dry and dust-free environment.

|                                                                |                      |
|----------------------------------------------------------------|----------------------|
| Minimum ambient temperature (without relative humidity)        | -25°C                |
| Maximum ambient temperature (with relative humidity up to 93%) | +70°C                |
| Humidity                                                       | max. 93%             |
| Atmospheric pressure                                           | 700 hPa to 1,060 hPa |

## Troubleshooting

### 1. Please refer to the list of solutions below when experiencing trouble with your Aerogen® Go nebuliser.

| Problem                                                             | Possible Cause                                                                                        | Solution                                                                                                                                                |
|---------------------------------------------------------------------|-------------------------------------------------------------------------------------------------------|---------------------------------------------------------------------------------------------------------------------------------------------------------|
| Weak or erratic aerosol stream.                                     | Medication reservoir is not properly attached.                                                        | Re-attach the medication reservoir as shown in pages 8-9 and power device on again.                                                                     |
|                                                                     | Medicine residue from previous inhalation sessions remains in the medication reservoir.               | Clean and disinfect the medication reservoir according to the procedures shown in page 10. When the device is fully cleaned and dry, power it on again. |
|                                                                     | Vibrating mesh component is clogged.                                                                  | Clean the medication reservoir. If the problem persists, please contact Aerogen Customer Service and replace the medication reservoir.                  |
|                                                                     | Electrodes on main unit or medication reservoir are wet or dirty.                                     | Clear the liquid or dirt and power the device on again.                                                                                                 |
|                                                                     | Electrodes on main unit or medication reservoir are blocked from touching.                            | Remove whatever is blocking the electrodes and power the device on again.                                                                               |
| The power indicator flashes briefly and then goes off.              | Medication reservoir is not properly attached.                                                        | Re-attach the medication reservoir as shown in pages 8-9 and power the device on again.                                                                 |
|                                                                     | No medicine in the medication reservoir.                                                              | Put liquid medicine in the medication reservoir.                                                                                                        |
|                                                                     | No contact between medicine and the mesh.                                                             | Adjust the nebuliser's angle relative to the user's mouth so the medicine contacts the vibrating mesh.                                                  |
|                                                                     | Electrodes on main unit and medication reservoir are wet or dirty.                                    | Clear the liquid or dirt and power the device on again.                                                                                                 |
| Power indicator is not lit and the device is not producing aerosol. | Batteries are installed backwards.                                                                    | Reinsert the batteries in the correct orientation and power the device on again.                                                                        |
|                                                                     | Connection between power adaptor and the device is incorrect.                                         | Remove and reconnect the power adaptor to the device by pushing the connector firmly into the electrical port on the device. Power on the device again. |
|                                                                     | Inner battery cover is not closed properly.                                                           | Twist off the outer battery cover and press down on the inner battery cover to seal it.                                                                 |
|                                                                     | Batteries are leaking.                                                                                | Carefully separate the medication reservoir from the main unit and discard the main unit.                                                               |
|                                                                     | Low battery indicator flashes; battery power is insufficient.                                         | Replace the batteries with new ones and power the device on again                                                                                       |
|                                                                     | Low battery indicator is lit constantly; battery power is insufficient; or battery power is depleted. | Replace the batteries with new ones and power the device on again.                                                                                      |

# Troubleshooting

| Problem                                                                           | Possible Cause                                                                                               | Solution                                                                                                                               |
|-----------------------------------------------------------------------------------|--------------------------------------------------------------------------------------------------------------|----------------------------------------------------------------------------------------------------------------------------------------|
| Power indicator is lit and the device is not producing aerosol.                   | Vibrating mesh is damaged or destroyed.                                                                      | Contact Aerogen Customer Service and replace the medication reservoir.                                                                 |
|                                                                                   | Electrodes on main unit or medication reservoir are wet or dirty.                                            | Clear the electrodes of liquid or dirt and power the device on again.                                                                  |
|                                                                                   | Vibrating mesh component is clogged.                                                                         | Clean the medication reservoir. If the problem persists, please contact Aerogen Customer Service and replace the medication reservoir. |
| No aerosol is generated while there is medicine left in the medication reservoir. | Medication reservoir is loose or otherwise not properly attached.                                            | Re-attach the medication reservoir as shown in pages 8-9 and power the device on again.                                                |
|                                                                                   | Connection between power adaptor and the device is incorrect.                                                | Remove and reconnect power adaptor from the device and power it on again.                                                              |
|                                                                                   | There are bubbles on the vibrating mesh component.                                                           | Shake the medication reservoir lightly to remove the bubbles and power the device on again.                                            |
|                                                                                   | The user changes the orientation of the device held in his or her hand while in use.                         | Hold the device vertically as the figure on page 9 and do not tilt it away from the user.                                              |
|                                                                                   | Medication reservoir is damaged or destroyed.                                                                | Contact Aerogen Customer Service and replace the medication reservoir.                                                                 |
| Device does not shut off automatically when medicine is depleted.                 | Some liquid medicines may produce foam in the medication reservoir.                                          | Switch off the device. Clean out the foam and power the device on again.                                                               |
|                                                                                   | Electrodes on main unit or medication reservoir are wet or dirty.                                            | Switch off the device. Clean the electrodes of liquid or dirt and power the device on again.                                           |
|                                                                                   | Medication reservoir is damaged or destroyed.                                                                | Switch off the device. Contact Aerogen Customer Service and replace the medication reservoir.                                          |
| Liquid is leaking from the medication reservoir.                                  | The medication reservoir cap is not sealed.                                                                  | Twist the medication reservoir cap counterclockwise to open it, then twist clockwise until it stops to reseal it.                      |
|                                                                                   | The medication reservoir cap is damaged.                                                                     | Contact Aerogen Customer Service and replace the medication reservoir.                                                                 |
|                                                                                   | Medication reservoir is damaged or destroyed, or the waterproof seals have aged and are no longer effective. | Contact Aerogen Customer Service and replace the medication reservoir.                                                                 |
|                                                                                   | Vibrating mesh is damaged or destroyed.                                                                      | Contact Aerogen Customer Service and replace the medication reservoir.                                                                 |

**2. If your Aerogen® Go nebuliser still does not function after following the solutions above, please contact Aerogen Customer Service and quote the Lot No. of medication reservoir or the unique device identification (UDI) of the nebuliser (illustrated on page 13) for assistance.**

## Disposal

The Aerogen® Go complies with the Directive 2012/19/EU of the European Parliament and of the Council on waste electrical and electronic equipment (WEEE). All of the components must not be disposed of with domestic waste. Specific parts such as batteries, power adaptor, and main unit, must be disposed of ecologically and in accordance with local regulations. For further information regarding the waste disposable service, contact Aerogen Customer Service. Correct handling of live parts and recycling would reduce the use of raw materials and contribute to the environment sustainability.

## Spare Parts and Accessories

To replace parts or accessories, contact Aerogen Customer Service and indicate the description: medication reservoir, mouthpiece, adult mask, paediatric mask, and power adaptor. For trouble shooting and warranty service, please quote the Lot No. of medication reservoir or the unique device identification (UDI) of the nebuliser.

Lot No. of medication reservoir

UDI of the nebuliser

## Product Specifications

|                     |                                                                                                                                                                                                          |
|---------------------|----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| Name                | Aerogen® Go Vibrating Mesh Nebulizer                                                                                                                                                                     |
| Model               | HCM-86C                                                                                                                                                                                                  |
| Dimensions          | 74 mm (L) × 46 mm (W) × 96 mm (H)                                                                                                                                                                        |
| Weight              | Approx. 75g (without batteries)                                                                                                                                                                          |
| Power supply        | 2×1.5 V AA alkaline batteries or the power adaptor.<br>The type C power adaptor (UMVUE3010-033020SA) or type G power adaptor (UMVUK3010-033020SA) has an input of 100-240 volts and output of 3.3 volts. |
| Power consumption   | Approx. 1.2W                                                                                                                                                                                             |
| Min. fill volume    | 1ml                                                                                                                                                                                                      |
| Max. fill volume    | 10ml                                                                                                                                                                                                     |
| Sound               | ≤ 50 dB                                                                                                                                                                                                  |
| Operating condition | 5°C~40°C, 15~93%R.H. (non-condensing), 700-1060hPa                                                                                                                                                       |

The Aerogen® Go Vibrating Mesh Nebulizer complies with the provisions of the Regulation (EU) 2017/745 (Medical Device Regulation) and the European Standard EN ISO 27427:2019.


## Technical Data

The particle size is measured according to EN ISO 27427:2019

Particle size: MMAD  $2.04 \pm 0.05 \mu\text{m}$

Aerosol output rate:  $0.7 \pm 0.07 \text{ ml/min}$  (2ml NaF 2.5%)

Aerosol output:  $1.54 \pm 0.15 \text{ ml}$  (2ml NaF 2.5%)

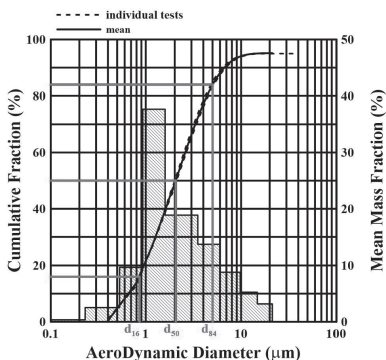The essential performance of the Aerogen® Go is illustrated as above. Performance may vary among different medicines such as suspension or high viscosity.

# Electromagnetic Compatibility

## Important information regarding Electro Magnetic Compatibility (EMC)

The Aerogen® Go conforms to EN60601-1-2:2015/A1:2021 standard for both immunity and emissions. However, certain devices such as mobile phones or other devices, which generate strong electrical or electromagnetic fields, may interrupt operation of the Aerogen® Go when they are used near the device.

### Manufacturer's Declaration – Electromagnetic Emissions

The Aerogen® Go is intended for use in the electromagnetic environment specified below. The customer or the user of the Aerogen® Go should assure that it is used in such environment.

| Emission Test                                                | Compliance | Electromagnetic Environment Guidance                                                                                                                                                                                    |
|--------------------------------------------------------------|------------|-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| RF Emissions<br>CISPR<br>11:2015+A1:2016+A2:2019             | Group 1    | The device uses RF energy only for its internal function. Therefore, its RF emission is very low and not likely to cause any interference nearby electronic equipment.                                                  |
| Harmonic Current Emissions<br>IEC 6100-3-2:2014              | Class B    | The device is suitable for use in all establishments, including domestic establishments and those directly connected to the public low voltage power supply network that supplies buildings used for domestic purposes. |
| Voltage Fluctuations & Flicker Emissions<br>EC 6100-3-3:2013 | Class A    |                                                                                                                                                                                                                         |
|                                                              | Complies   |                                                                                                                                                                                                                         |

### Manufacturer's Declaration – Electromagnetic Immunity

The Aerogen® Go is intended for use in the electromagnetic environment specified below. The customer or the user of the Aerogen® Go should assure that it is used in such environment.

| Immunity Test                                         | IEC 60601 Test Level                                                                                                                                                                    | Compliance | Electromagnetic Environment Guidance                                                                                                                                                                                                                                                                                                |
|-------------------------------------------------------|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|------------|-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| Electrostatic Discharge (ESD)<br>IEC 61000-4-2:2008   | Contact: $\pm 8\text{kV}$<br>Air: $\pm 2\text{kV}$ , $\pm 4\text{kV}$ ,<br>$\pm 8\text{kV}$ , $\pm 15\text{kV}$                                                                         | Complies   | Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30%.                                                                                                                                                                                       |
| Electrical Fast Transient/Burst<br>IEC 61000-4-4:2012 | $\pm 2\text{kV}$ for AC and<br>DC 5/50 Tr/Th<br>ns100kHz                                                                                                                                | Complies   | Mains power quality should be that a typical professional healthcare facility environment or home healthcare environment.                                                                                                                                                                                                           |
| Surge<br>IEC 61000-4-5:2014+A1:2017                   | $\pm 0.5\text{kV}$ , $\pm 1.0\text{kV}$<br>line(s) to line(s)<br>$\pm 0.5\text{kV}$ , $\pm 1\text{kV}$ ,<br>$\pm 2\text{V}$ line(s) to<br>earth<br>1.2/50 (8/20)<br>Tr/Th $\mu\text{s}$ | Complies   | Mains power quality should be that a typical professional healthcare facility environment or home healthcare environment.                                                                                                                                                                                                           |
| Power Frequency Magnetic Field<br>IEC 61000-4-8:2009  | 30 A/m, 50Hz or<br>60Hz                                                                                                                                                                 | Complies   | Power frequency magnetic fields should be at levels characteristic of a typical location in a typical professional healthcare facility environment or home healthcare environment.<br>Keep the device away from sources of high level of power line magnetic fields (in excess of 30 A/m) to reduce the likelihood of interference. |

# Electromagnetic Compatibility

| Immunity Test                                                                                 | IEC 60601 Test Level                                                                                                                                  | Compliance | Electromagnetic Environment Guidance                                                                                                                                                                                                                                                                                                                                                                                                                                                                  |
|-----------------------------------------------------------------------------------------------|-------------------------------------------------------------------------------------------------------------------------------------------------------|------------|-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| Voltage Dips, Short Interruptions, and Voltage Variations<br><br>IEC 61000-4-11: 2004+A1:2017 | Voltage dips:<br>0% UT for 0.5 cycle<br>0% UT for 1 cycle<br>70% UT for 25/30 cycles (50/60 Hz)<br>Voltage interruptions:<br>0% UT for 250/300 cycles | Complies   | Mains power quality should be that a typical professional healthcare facility environment or home healthcare environment. If the user requires continued operation during power mains interruptions, it is recommended that the device be powered from an uninterruptible power supply or batteries.                                                                                                                                                                                                  |
| Conducted Disturbances Induced by Radio Frequency Fields<br><br>IEC 61000-4-6:2013            | 3Vrms:<br>150 kHz-80 MHz<br>6Vrms: In ISM bands between 150 kHz and 80 MHz<br>80% AM at 1 kHz                                                         | Complies   | Portable and mobile RF communications equipment should be used no closer to any part of the unit, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter as below.                                                                                                                                                                                                                                                    |
| Radiated, Radio Frequency, Electromagnetic Field<br><br>IEC 61000-4-3:2010                    | 10 V/m<br>80 MHz-2.7 GHz<br>80% AM at 1 kHz                                                                                                           | Complies   | <b>Recommended separation distance</b><br>$d=1.17\sqrt{P}$<br>$d=0.35\sqrt{P}$ 80 MHz to 800 MHz<br>$d=0.7\sqrt{P}$ 800 MHz to 2.7 GHz<br>where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in meters (m).<br>Interference may occur in the vicinity of equipment market with the following symbol:  |

**Note 1:** At 80 MHz and 800 MHz, the higher frequency range applies.

**Note 2:** These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects, and people.

## Recommended separation distances between portable and mobile RF communications equipment and the Aerogen® Go.

The Aerogen® Go is intended for use in the electromagnetic environment in which radiated RF disturbance are controlled. The customer or the user of the Aerogen® Go can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communication equipment (transmitters) and the Aerogen® Go as recommended below, according to the maximum output power of the communication equipment.

| Rated maximum output power of transmitter W | Separation distance according to frequency of transmitter M |                                      |                                       |
|---------------------------------------------|-------------------------------------------------------------|--------------------------------------|---------------------------------------|
|                                             | 150 kHz to 80 MHz<br>$d=1.2\sqrt{P}$                        | 80 MHz to 800 MHz<br>$d=1.2\sqrt{P}$ | 800 MHz to 2.5 GHz<br>$d=2.3\sqrt{P}$ |
| 0.01                                        | 0.12                                                        | 0.12                                 | 0.23                                  |
| 0.1                                         | 0.38                                                        | 0.38                                 | 0.73                                  |
| 1                                           | 1.2                                                         | 1.2                                  | 2.3                                   |
| 10                                          | 3.8                                                         | 3.8                                  | 7.3                                   |
| 100                                         | 12                                                          | 12                                   | 23                                    |

For transmitters rated at a maximum output power not listed above, the recommended separation distance d in meters (m) can be estimated using the equation applicable to the frequency of the transmitter, where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.

**Note 1:** At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies.

**Note 2:** These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects, and people.

# Symbols and Machine Markings

## Symbols

|                                                                                   |                                                       |                                                                                   |                                                                  |
|-----------------------------------------------------------------------------------|-------------------------------------------------------|-----------------------------------------------------------------------------------|------------------------------------------------------------------|
|  | Follow instructions for use                           |  | Catalogue number                                                 |
|  | Caution                                               |  | Lot number                                                       |
|  | Manufacturer                                          |  | The device is protected against water jets and dust              |
|  | Date of manufacture                                   |  | Minimum and maximum ambient temperature                          |
|  | Distributor                                           |  | Minimum and maximum humidity                                     |
|  | Importer                                              |  | Minimum and maximum air pressure                                 |
|  | This product complies with 93/42/EEC (Medical device) |  | Disposal                                                         |
|  | Authorised representative in the European community   |  | Protection class of the part used: Type BF (Masks or mouthpiece) |
|  | Authorised representative in Switzerland              |  | Single patient multiple use                                      |
|  | This product complies with the UK MDR 2002            |  | Keep dry                                                         |
|  | Unique device identifier                              |  | Fragile                                                          |
|  | Quantity                                              |  | Consult instructions for use                                     |

## Machine markings

|                                                                                     |                 |                                                                                     |                       |
|-------------------------------------------------------------------------------------|-----------------|-------------------------------------------------------------------------------------|-----------------------|
|  | On/Off button   |  | Low battery indicator |
|  | Power indicator |                                                                                     |                       |

# Warranty Information

## Warranty Information – Limited Warranties

The Aerogen® Go Vibrating Mesh Nebulizer, excluding the medication reservoir, is warranted to be free from defects in materials and workmanship within one (1) year from the date of purchase when used in accordance with the instruction manual provided with the device.

The warranty extends only to the original purchaser. Aerogen will, within the reasonable options, replace without charge your Vibrating Mesh Nebulizer. Replacement is Aerogen's only responsibility and your only remedy under the above warranty.

To obtain warranty service, please contact Aerogen Customer Service at [info@aerogen.com](mailto:info@aerogen.com) or +353 91 540 400.

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Aerogen Ltd.
HCmed Innovations Co., Ltd.
Medical Device Safety Service GmbH.,
Medical Device Safety Service CH GmbH.,
MDSS UKRP Ltd., 6 Wilmslow Road,